


THE TRAVelite

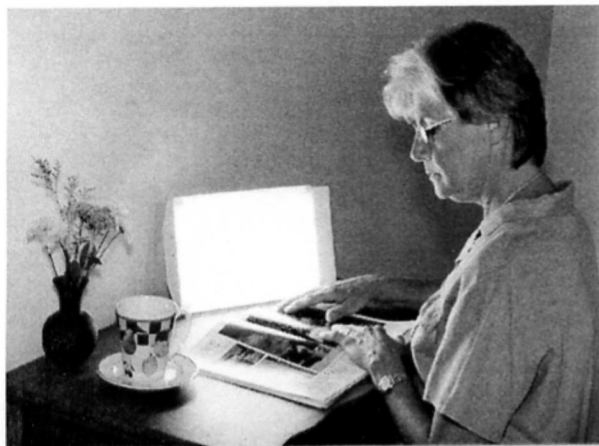

The TRAVelite is a compact and versatile unit that will let you enjoy the light in a variety of situations. It is a great travel light – it fits in a briefcase and weighs only 1 Kg (2 Lbs). The unit comes with a bracket that allows you to stand the TRAVelite vertically or horizontally. See back for more information. The unit can also be wall mounted or turned into a floor lamp by mounting it on a tripod. Now you can decide when the sun shines!

The TRAVelite delivers 6,000 Lux of light at a distance of 15 inches (40 cm.). We recommend that you use the lamp for about one hour each day, preferable first thing in the morning, sitting as shown in our brochure photo. You do not need to look directly at the light. The light needs to shine in your eyes. Do

not wear dark glasses while using light. Use the light regularly at about the same time each day. Also try to spend time outdoors during sunny periods. Eat sensibly, exercise, and try to keep regular wake – sleep hours.

PROBLEMS:

1) UNIT DOES NOT LIGHT:

If the unit does not light 3 seconds after the switch is turned ON, please check that the plug is well inserted into the outlet. If this fails, try other outlets. If the problem persists, remove the plastic diffuser by gently removing the four nylon screws. Check that both fluorescent tubes are equally well seated in the sockets. Also, if there is an area of the tube that is clear glass, without the white powder coating on the inside, the tube is defective. If your lamp has just been brought in from the cold, you must wait until it reaches room temperature before using it. Still no luck, call us. We will try to help you.

2) OTHER PROBLEMS:

Light can cause eye and/or skin irritation and headaches. This is most prevalent in people with blue or green eyes, blonde or red hair, and/or sensitive skin. All these problems can usually be resolved by gradually increasing use of lamp. First, stop using light until eye or skin irritation clears up. Then, start using light five minutes a day. Maybe increase the distance between you and the lamp as well. Gradually, over a period of about two weeks, you should be able to lengthen your exposure and shorten the distance and comply with your original light regimen. If you still have a problem, call us, we may be able to help you.

Please note that some medications such as lithium, melatonin, tetracycline, St. John's Wort, phenothiazines, hematoporphyrins, acne creams, creams with retinoic acid, chloroquine, and other chemotherapy drugs can make you very sensitive to light.

MAINTENANCE:

The lamp should not require maintenance other than occasional dusting. Smudges can be wiped with a soft damp cloth. The fluorescent tubes emit less light as they get older. Typically, after 2,000 hours of usage they emit 15% less light.

WARRANTY:

We undertake to repair or replace a defective lamp at no charge for a period of five years. Shipping costs to our plant must be paid by you. You **MUST** call us at **800 263-0066** before returning a unit.



GUARANTEE:

Not all vendors have the same merchandise return policy. Only the conditions at the point of purchase apply.

WARNING:

If you have a health problem, if you have been told you must wear sunglasses in bright light, or if you have a history of eye problems, please consult a health-care professional before using our lamp.

Light intensity of the TRAVelite	Distance from the lamp	Suggested usage per day with the TRAVelite
6000 lux	15 inches (40 cm)	60 minutes
10,000 lux	10 inches (25 cm)	30 minutes

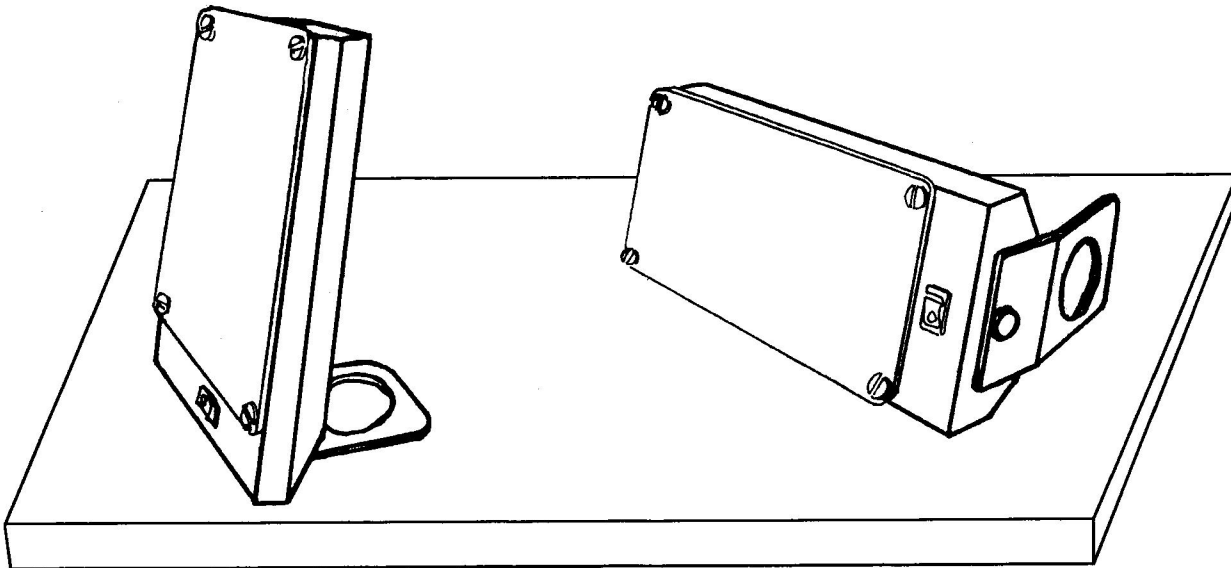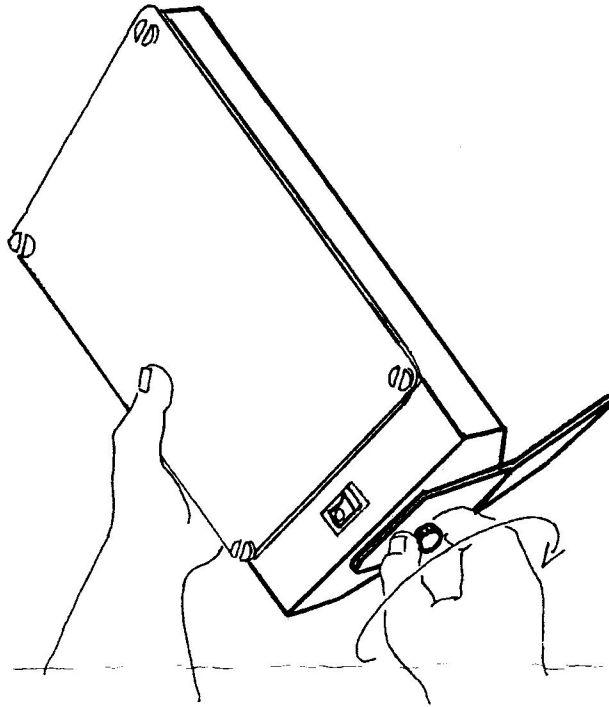

Northern Light Technologies 8971 boul. Henri-Bourassa West, Montreal (Quebec) H4S 1P7 Canada

Tel: (514) 335-1763 Toll free: 1 800 263-0066 Fax: (514) 335-7764

Website: www.northernlight-tech.com Email: info@northernlight-tech.com